

Dr.Yglo

WART & PLANTAR WART REMOVER

INDICATIONS FOR USE: DR. YGLO® WART REMOVER IS INDICATED FOR OVER-THE-COUNTER TREATMENT OF COMMON WARTS AND PLANTAR WARTS FOR PATIENTS AGE 4 YEARS AND OLDER

PLEASE READ INSERT CAREFULLY BEFORE USING DR. YGLO® WART AND VERRUCA TREATMENT. WE ADVISE YOU TO KEEP THIS LEAFLET, YOU MAY NEED TO READ IT AGAIN.

TABLE OF CONTENTS

1. CONTENTS OF DR.YGLO WART REMOVER
2. BACKGROUND ON WARTS
3. WHAT IS DR.YGLO WART REMOVER USED FOR?
4. WARNINGS
5. CAUTIONS
6. WHAT TO DO IF YOU ARE UNCERTAIN IF THE CONDITION IS A WART?
7. HOW DOES DR.YGLO WART REMOVER WORK?
8. SIDE EFFECTS
9. INSTRUCTIONS FOR USE
10. APPLICATION TIMES
11. FOLLOW-UP TREATMENT

1. CONTENTS OF DR.YGLO WART REMOVER

A pressurized canister filled with 50ml of the compressed liquid gas dimethyl ether. A custom application unit used to administer the cold delivered by the cryogen to the wart. This information leaflet in which you can read about the product how it works, how to use the product to treat warts, warnings and limitations. A planter Wart or Verrucas are similar to warts but occur only on the sole of the foot or underneath the toes. Generally painful, they often look like a small white ring of skin with black dots in the centre. Verrucas are also contagious and can be transmitted to another person by contact with a wet surface such as a floor.

2. BACKGROUND ON WARTS

Warts are growths on the skin or mucous membranes (such as inside mouth, nose, anus, genitals, lips) that are caused by the human papilloma virus (HPV). There are more than 80 different types of HPV. Many are attracted to the skin where they cause different skin-related warts including plantar warts and common warts. Others are attracted to the mucus membranes where they can cause warts on the genitals or anus. Skin-related warts such as common warts and plantar warts are generally harmless. However, genital warts have been linked to cancer in the cervix. Warts affect 7 to 10% of the general population. They are most common in children and adolescents. Warts are very contagious and may appear more often if the body's defense system is weak. Warts are passed from person to person or from one body location to another on the same person by touch or contact with flakes of skin. More than 50% of warts disappear within two years with no treatment. However, they are unsightly and can cause pain in the case of plantar warts. Therefore, patients often choose to treat warts rather than wait for them to disappear. Treatments for warts include different topical acid preparations, surgical removal, burning or cryotherapy (freezing). Duct tape has also recently been reported in the news to be effective. Liquid nitrogen is the most commonly applied freezing method, but it is only available through the physician. Dr. Yglo Wart Remover is an alternative to liquid nitrogen. It contains dimethyl ether and is as effective and well tolerated as liquid nitrogen.

3. WHAT IS DR.YGLO WART REMOVER USED FOR?

Dr. Yglo Wart Remover is used for the removal of common warts as well as the removal of plantar warts on the bottom of the foot. Common warts and plantar warts are harmless growths that consist of a center of inner skin (or connective tissue) with blood vessels covered with layers of outer skin (or epidermis tissue). Common warts are skin-coloured, usually appear on the hands and feet and are easily recognized by the rough "cauliflower-like" appearance of the surface. Plantar warts only occur on the bottom of the foot and are often flat due to the pressure caused by walking. Plantar warts grow inward, are often calloused, and can be very painful. The footprint pattern is smooth over the area where the plantar wart occurs because "skin lines" do not grow across the surface of plantar warts.

4. WARNINGS

- Do not use if you cannot see clearly to read the information on the container or this Information Leaflet. Use only as directed.
- Do not use if you cannot clearly see the wart you are treating. If you mistakenly apply Dr. Yglo Wart Remover on the conditions that are not warts, it may cause serious burns and permanent scarring of the skin.
- Do not use on children under the age of 4 years. It is uncertain how Dr. Yglo Wart Remover may affect the child's skin.
- Do not use if you are diabetic or if you have poor blood circulation. Dr. Yglo Wart Remover works by freezing the wart. The frozen area in diabetic patient and those with poor blood circulation may not heal as quick or completely.
- Do not use if you are pregnant or breast-feeding. It is uncertain how Dr. Yglo Wart Remover may affect patients under these conditions.
- Do not use on area with thin skin such, face, armpits, breasts, bottom, or genital area. Use on these areas may cause burns and permanent scarring.
- Do not use on the irritated skin, or on any area that is infected, reddened or showing any signs of inflammation, such as itching or swelling. Treating these areas with Dr. Yglo Wart Remover will result in worsening of the condition and possible further infection. These conditions may not be wart and should be examined by doctor.
- Do not use on moles, birthmarks, warts with hair growing from them, genital (on penis or vaginal area) warts, or warts on the face or mucus membranes (such as inside mouth, nose, anus, genital, lips, ears or near eyes). Use on these areas may cause burns and permanent scarring of the skin or blindness.

Use Dr. Yglo Wart Remover only if you are sure the skin condition is a wart.

If you are not sure, consult your doctor. If you treat a skin condition as a wart with Dr. Yglo Wart Remover and it is not a wart, you delay proper treatment. This delay may lead to a more serious skin condition or disease.

Avoid Contact with Eyes. Contact with eyes may cause blindness. Do not rub or touch eyes with hands if Dr. Yglo Wart Remover has been applied to them. If product gets into eyes, flush with water for 15 minutes and call medical help right away.

Do not Swallow. Content is toxic and may cause serious internal damage. If swallowed, get medical help or contact a Poison Control Center right away.


For External Use Only. Use only as directed. If Dr. Yglo Wart Remover is not used exactly as instructed or you mistakenly apply it directly to the skin or use it on conditions that are not warts, it may cause serious burns and permanent scarring skin.

Keep out or Reach of Children. Children cannot understand the instructions to properly use this product. Content is hazardous if swallowed, inhaled, sprayed on the parts other than wart, or exposed to heat.

5. CAUTIONS

- The lot number and the expiration date appear on the bottom of the aerosol can as well as the outside of the carton.
- Do not use Dr. Yglo Wart Remover if the cone is stuck to position "1".
- If being used by more than one person and/or after each application, clean the applicator with alcohol on a small piece of cloth.
- If there is more than one wart on a finger or toe, only treat one wart at a time. Freezing more than one wart on a single finger or toe is too large of an area to freeze at one time and may take longer to heal or the area may not heal properly. Treat each remaining wart on a single finger or toe separately, with a two-week period between treatments.
- In most cases, only one treatment is needed for each wart. Warts are usually gone in two weeks after treatment with Dr. Yglo Wart Remover. If a wart or a part of it is still there two weeks after treatment, you may then safely treat again. Do not treat each wart with Dr. Yglo Wart Remover more than three times in total. Consult your doctor if three treatments have not led to improvement.
- The aching/stinging feeling in treated common warts should disappear

HOW TO USE DR YGLO WART REMOVER

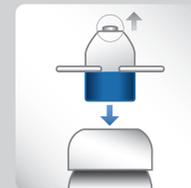


1. Apply the actuator on the aerosol can, with the metal tip facing up. Do not push down on the activator until ready to cool.



2. Place the aerosol can on a hard surface and push down on the white actuator for 3 seconds. Then allow the can to rest for 2-3 seconds. Be careful to avoid contact with the metal tip when pressing down.



3. Remove the actuator from the aerosol can, with the metal tip facing up.

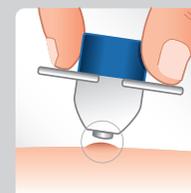

4. Apply the metal tip on the wart for 20 seconds or verruca for 40 seconds. Make sure the metal tip covers as much of the wart or verruca as possible.

Dr.Yglo

WART & PLANTAR WART REMOVER

within few hours after application. For plantar warts, the aching/ stinging feeling may last 24 hours. If discomfort or aching/ stinging feeling beyond this time, 3 hours for common warts and 24 hours for plantar warts, talk to your doctor. This may mean that the area has been frozen too deeply and medical attention may be needed.

- Do not inhale vapour/spray and use only in well-ventilated areas because it is uncertain how the vapour/spray will affect the lungs. Breathing vapour/ spray may be toxic.
- Wash hands (except treated area) immediately after use. Washing the treated area immediately may be uncomfortable since you may still be experiencing some aching/stinging. The treated area may be washed after the aching/stinging subsides.

Extremely Flammable. Content under pressure. Do not pierce or burn aerosol spray can, even after use. Protect from sunlight and do not expose to temperature above 120°F (50°C). Store at room temperature away from heat. Keep away from fire or flame and from source of sparks or ignition. Do not smoke while using Dr. Yglo Wart Remover. Piercing the container or exposure to the heat may cause the container to explode, causing serious injury. Do not breathe the vapour/spray and use only in well-ventilated areas because it is uncertain how the vapour/spray will affect the lungs. Breathing vapour/spray may be toxic.

6. WHAT TO DO IF YOU ARE UNCERTAIN IF THE CONDITION IS A WART?

See "Warnings" and "Cautions." Use Dr. Yglo Wart Remover only if you are sure the skin condition is a wart. If you are not sure, consult your doctor. If you treat a skin condition as a wart with Dr. Yglo Wart Remover and it is not a wart, you delay proper treatment. This delay may lead to a more serious skin condition or disease.

7. HOW DOES DR YGLO WART REMOVER WORK?

1. Dr. Yglo Wart Remover freezes warts on the spot. The freezing kills both the virus and the infected skin. When the applicator is pressed (applied) on the wart, freezing begins within a few seconds and the skin will turn white (see diagram and instructions for use). A slight aching, itching and/or stinging feeling will occur as a result of the freezing. After the applicator is removed, colour will return and the aching, itching, stinging feeling will quickly decrease, disappearing within a few hours.
2. After the area thaws and colour returns, a redness that is the same size as the frozen region may appear. As a result of freezing, a blister may form (sometimes filled with blood) under the wart within a few days. You may not see this easily.
3. The frozen skin with the treated wart will fall off after about 10-14 days. In the meantime, new skin will have formed under the blister.

8. SIDE EFFECTS

- A stinging, itching, aching feeling during and after freezing that quickly decreases after the area thaws.
- Possible changes in colouring from white to red.
- Possible blistering.
- Black dots may appear in the middle of the warts after a few days.

9. INSTRUCTIONS FOR USE (COMMON WARTS)

1. Apply the actuator on the aerosol can, with the metal tip facing up. Do not push down on the activator until ready to cool.
2. Place the aerosol can on a hard surface and push down on the white actuator for 3 seconds. Then allow the can to rest for 2-3 seconds. Be careful to avoid contact with the metal tip when pressing down.
3. Remove the actuator from the aerosol can, with the metal tip facing up.
4. Apply the metal tip on the wart for 20 seconds. Make sure the metal tip covers as much of the wart or verruca as possible.

10. INSTRUCTIONS FOR USE (COMMON WARTS)

Soak foot with plantar wart for 5 min. in warm water. Dry thoroughly. File the plantar wart area with a pumice stone or pumice file (not included in the Dr. Yglo Wart Remover (packaging) before applying Dr. Yglo Wart Remover. Follow the instruction as indicated at point 1,2,3 below and activate Dr. Yglo Wart Remover.

1. Apply the actuator on the aerosol can, with the metal tip facing up. Do not push down on the activator until ready to cool.

2. Place the aerosol can on a hard surface and push down on the white actuator for 3 seconds. Then allow the can to rest for 2-3 seconds. Be careful to avoid contact with the metal tip when pressing down.
3. Remove the actuator from the aerosol can, with the metal tip facing up.
4. Apply the metal tip on the wart for verruca for 40 seconds. Make sure the metal tip covers as much of the wart or verruca as possible.

Comfort pads (not included in the Dr. Yglo Wart Remover packaging) can be used to protect the tender plantar wart area after treatment as been completed. Carefully remove pad from backing. Apply adhesive side of pad to skin with hole directly over the small tender spot.

11. INSTRUCTIONS FOR USE (COMMON AND PLANTAR WARTS)

For each common wart, only one treatment is necessary in most cases. Plantar warts may require more than one treatment. In two weeks after treatment with Dr. Yglo Wart Remover, common warts are usually gone. If a Wart or part of it is still there two weeks after treatment, you may then safely treat it again. Do not treat each wart with Dr. Yglo Wart Remover more than three times in total. Consult your doctor if three treatments have not led to improvement. Do not apply Dr. Yglo Wart Remover together with any other method of wart removal because it is uncertain how the combination may affect your skin. Combination usage may result in serious burns and permanent scarring of the skin.

It is important that the instructions in the Information Leaflet are followed exactly. Careless use can cause freezing too deeply, causing permanent damage to the uppermost layer of the skin on fingers and toes. This may result in permanent scars or nerve damage.

12. INSTRUCTIONS FOR USE (COMMON AND PLANTAR WARTS)

Precisely apply the frozen metal tip on the common wart (hand) for 20 seconds or less, depending on the size of the wart. If you use it on a plantar wart on an area of the foot with thicker or calloused skin, such as the balls of the foot or the heel, hold the applicator on the plantar wart for 40 seconds or less. Treat plantar warts occurring on less calloused areas of the foot, such as the arch of the foot or toes, for 20 seconds or less. Calloused areas of the foot vary from person to person. If the area of the foot where the plantar wart occurs is not calloused, do not apply Dr. Yglo Wart Remover for more than 20 seconds.

12. FOLLOW-UP TREATMENT

- After application, the colour of the skin may change, varying from white to red. A blister may form under the wart and an aching/stinging feeling normally will occur for up to a few hours. If the aching, stinging feeling continues until the next day, or if another problem (such as freezing the skin beyond the size of the wart or less of feeling) occurs, contact your doctor. It is normal for the treated area to be sensitive for a few days.
- Black dots may appear in the center of the treated wart after a few days.
- Keep the treated area clean.
- You may swim or shower.
- Do not pick or scratch the treated area as this may result in infection.
- Protect blisters, if necessary, with a gauze or sterile adhesive bandage.
- Do not puncture blisters. Puncturing blisters will be painful and may result in infection.

Dr. Yglo Wart Remover contains:


- Dimethyl ether

Last revision: January 2021

For more information and an instructional video please visit www.dryglo.com

Manufacturer / Distributor:
TheOTCLab
Fred. Roeskestraat 115
1076 EE Amsterdam
The Netherlands


HOW TO USE DR YGLO WART REMOVER

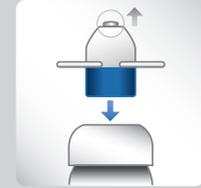

1. Apply the actuator on the aerosol can, with the metal tip facing up. Do not push down on the activator until ready to cool.



2. Place the aerosol can on a hard surface and push down on the white actuator for 3 seconds. Then allow the can to rest for 2-3 seconds. Be careful to avoid contact with the metal tip when pressing down.

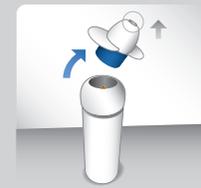

3. Remove the actuator from the aerosol can, with the metal tip facing up.



4. Apply the metal tip on the wart for 20 seconds or verruca for 40 seconds. Make sure the metal tip covers as much of the wart or verruca as possible.